

# Motif Luna® Quick Guide


## 1. Buttons & Features

### LCD Screen

- A Massage Mode
- B Expression Mode
- C Cycles
- D Vacuum Level
- E Timer

### Button Control


- F Mode Button
- G Cycles Increase/Decrease Button
- H LED Night Light
- I Levels Increase/Decrease Button
- J Power Button

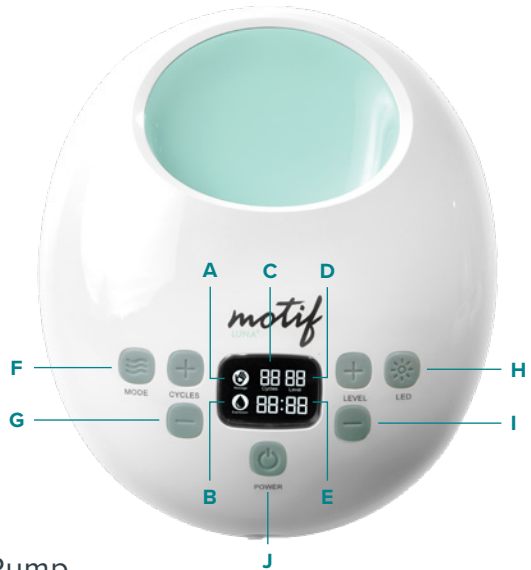

## 2. Preparing The Luna Pump


### Assembling The Luna Pump

1. Wash hands before touching breast pump components.
2. Insert the diaphragm into the diaphragm bottom (the component that connects to the breast shield). Ensure that it fits securely and that the diaphragm edge fits around the rim of the diaphragm bottom.
3. Connect the diaphragm top (the component that connects to the tubing) over the diaphragm bottom, and diaphragm, by snapping on the diaphragm bottom into the rim of the top diaphragm.
4. Push the valve onto the underside of the breast shield.
5. Connect the diaphragm assembly to the top of the breast shield. Use a twisting motion to attach, and later remove, this assembly to avoid damaging it.
6. Screw the bottle onto the underside of the breast shield.
7. Connect the tubing to the diaphragm assembly, and connect the other end of the tubing to the Luna. (For single-pumping attach the cap on one of the air flow ports).
8. Attach the power adapter to the Luna.
9. Plug the power adapter into the electrical outlet.

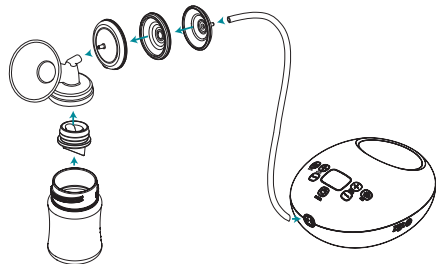


### 3. Using The Luna

#### HOW TO USE THE DUAL PUMP

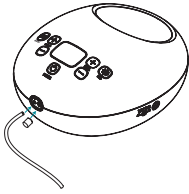

##### Single Pump

When changing back to Single-Pump mode make sure to cover the air flow port with the plug.

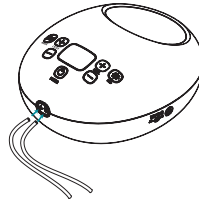

##### Double Pump

To use the Luna in Double-Pump mode remove the plug from the air flow port and connect the tubing from the extra collection kit.

#### TURNING ON/OFF THE LUNA PUMP

1. **Press the Power button.** The device will power ON, indicated by the screen lighting up. The Luna will begin in Massage Mode to help with faster letdown.
2. **Press the Mode button** to change to Expression Mode once your milk starts flowing. Massage Mode also has a function in which it will automatically end and switch to Expression Mode without having to press the button to change the modes manually.
3. **Adjust suction to preferred comfort level.**
4. **When you wish to stop expressing,** first press the Power button, then remove the breast shield and milk collection kit from your breast.

### 4. Warnings & Precautions



**WARNING:** Potentially dangerous situation that could result in serious injury.



**CAUTION:** A possibly dangerous situation that could result in slight bodily injury.



Manufactured for:  
**Motif Medical®**, LLC  
8 The Green, Suite 5601  
Dover, DE 19901

Phone: (844) 272-8390  
Email: [info@motifmedical.com](mailto:info@motifmedical.com)  
Web: [www.motifmedical.com](http://www.motifmedical.com)

*motif*  
MEDICAL